



# Alliance

TECH MEDICAL

## The Uniquely NEW All Flow Trigger Nebulizer Bottle

The AllFLOW Trigger Nebulizer Bottle is an innovative two-in-one nebulizer that produces aerosol at the user's discretion without taking the continuing-mode aerosol mechanism away from the user. AT-NEB helps the prescribed medication to be delivered deep into patient's lungs and stays as the patient paced it: most important of all, without exposing caregivers and the environment to risk of cross contamination.



### AllFLOW™ Trigger Nebulizer Bottle

- Ergonomic Bottle Design
- Built in Expiration Filter
- Lock On – Continuing Mode
- Lock Off – Discretionary Operation
- Leak Free
- Breathe Enhanced Technology

Connect with us online



[www.alliancetechnical.com](http://www.alliancetechnical.com)

## Product Specifications:

Features	AT-NEB
Delivery Method	Triggered or Continuing Mode
Nebulizing Method	Triggered or Continuing Mode
Fill Lines	0.5ml, 1ml, 2ml, 3ml, 4ml, 5ml, 6ml
Operating Flow Rate	7-8 lpm, 50 psi
Air Containment	Yes, By the Trigger
Latex Free	Yes
Removable Mouthpiece	Yes
Operating Angle	Up to 45 degrees
Aerosol Performance	
Fine Particle Fraction (% < 4.8 microns)	84%
Mass Median Aerodynamic Diameter	2.5um
Environmental Aerosol Emission	4%

Test Equipment: Spraytec Laser Beam • Test Method: Breath Enhanced Mode Nebulization

## Ordering Information:

ATM PART#	Quantity	Description
4730050	Box of 50	Triggered bottle + Mouth piece + Oxygen Tubing
4730010	Box of 10	
4740050	Box of 50	Triggered bottle + Mouth piece + Oxygen Tubing+Filter
4740010	Box of 10	
4750150	Box of 50	Triggered bottle + Oxygen Tubing + Aerosol Adult Mask
4750110	Box of 10	
4760150	Box of 50	Triggered bottle + Oxygen Tubing + Aerosol Child Mask
4760110	Box of 10	


ATM PART#	Description
4700010	Filter (10pcs/pack)
4700250	Aerosol adult mask + L-connector
4700350	Aerosol child mask + L-connector
4700450	Vent-free adult mask + L-connector
4700550	Vent-free child mask + L-connector
4700000	Oxygen tubing




Ergonomic Bottle Design



Leak Free

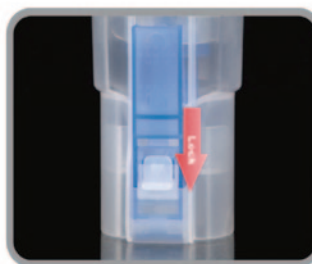

Lock On – Continuing Mode  
Lock Off – Discretionary Operation

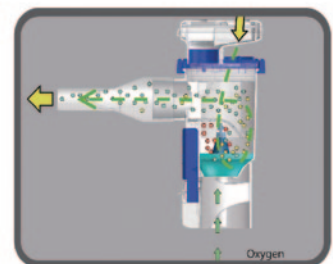

Breath Enhanced Technology (One Way Valve)



**Alliance Tech Medical, Inc.**  
 PO Box 6024 / Granbury, TX 76049  
 1.800.848.8923 / 817.326.6357  
 Fax: 817.326.2182 / [www.alliancetechnical.com](http://www.alliancetechnical.com)

**Eastern Office** 302.838.2638    **Central Office** 800.848.8923  
**Southern Office** 817.326.3183    **Western Office** 800.848.8923