




Peak Flow Meters

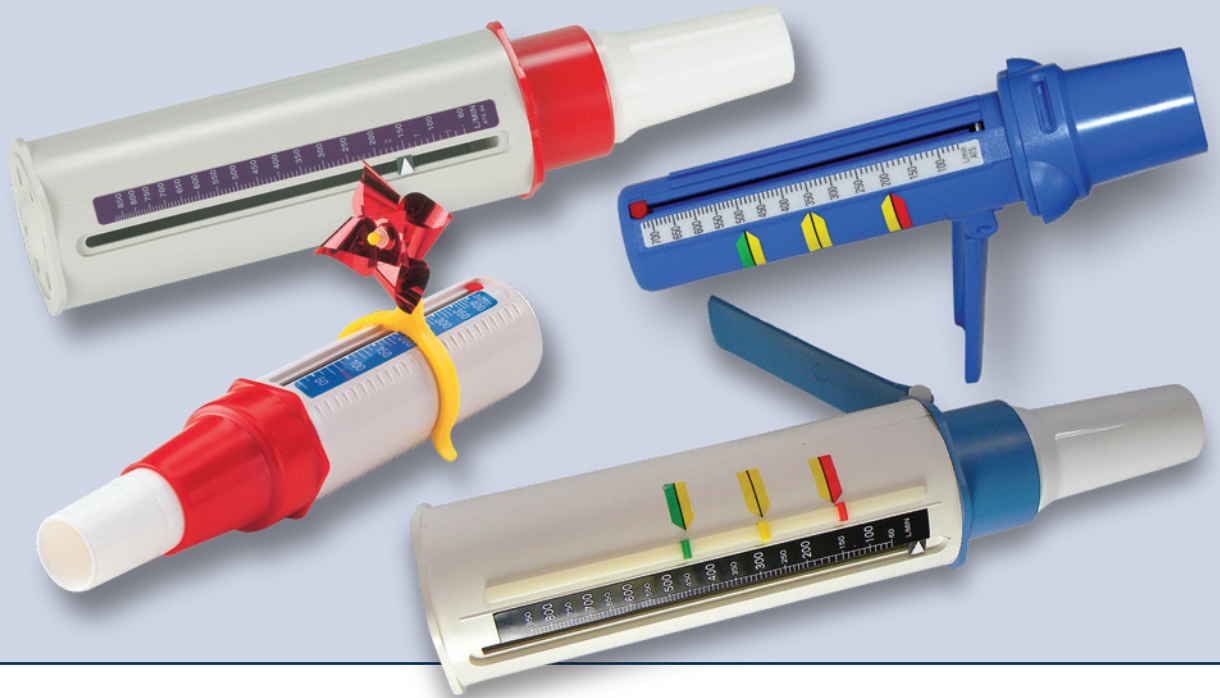

Patients can be confident in the proven, time tested performance of the Mini-Wright™ and Mini-Wright AFS™ peak flow meters, or the economical, accurate AllFLOW™ and Airzone™ peak flow meters.

Years of clinical tests have proven that both the Standard Range Mini-Wright™ and Low Range Mini-Wright AFS™ peak flow meters have delivered exceptional and consistent accuracy.

Expect the same from the AllFLOW™ and Airzone™ peak flow meters, which provide patients an economical choice of accurate peak flow monitoring.

- All meters for use by children and adults
- All individually calibrated
- Mini-Wright and AllFLOW models have built-in, one-way valves
- Optional Windmill Trainers for Mini-Wright models, teach proper peak flow technique and improve monitoring compliance. Ideal for children as young as 2-1/2 years
- AllFLOW and Airzone's folding handles adds to portability and convenience
- AllFLOW and Airzone's tape zone markers prevent inadvertent repositioning

Choose accurate peak flow monitoring as part of your asthma management program.

Manage your asthma



AllFLOW™ Specifications

ATS Scale	60-880 L/min
Scale Increments	10LPM
Size	5-7/8"H X 1-7/8"W
Weight	2.9 oz
Materials	All Latex Free
Patient Information	Operating Instructions/ Patient Diary in English, Spanish and French



Airzone™ Specifications

ATS Scale	60-720 L/min
Scale Increments	10LPM
Size	5-1/2"H X 1-3/4"W
Weight	1.5 oz
Materials	All Latex Free
Patient Information	Operating Instructions/ Patient Diary in English, Spanish and French



Mini-Wright™ and Mini-Wright AFS™ Specifications

ATS Scale	
MINI-WRIGHT Standard	60-880 L/min
MINI-WRIGHT AFS Low	30-400 L/min
Scale Increments	10LPM
Size	
MINI-WRIGHT Standard	5-7/8"H X 1-7/8"W
MINI-WRIGHT AFS Low	5-7/8"H X 1-1/2"W

Weight	
MINI-WRIGHT Standard	2.9 oz
MINI-WRIGHT AFS Low	2.2 oz

Materials All Latex Free

Patient Information Operating Instructions/
Patient Diary in English,
Spanish and French



PEAK FLOW METER

Disposable Mouthpieces

Economical and convenient. Reduce risk of cross-infection, helping to protect patient and equipment of contamination. One-way valve mouthpieces available for expiratory use.



Windmill Trainer

The Windmill Trainer is a revolution in teaching proper peak flow technique. By improving technique, use of the Windmill Trainer reduces variability in peak flow readings and improves monitoring compliance. Ideal for use with children as young as 2-1/2 years old, the elderly and those with language or learning difficulties. The Windmill Trainer is an effective yet enjoyable tool. Use the Windmill Trainer with both the Mini-Wright AFS Low Range and Mini-Wright Standard Range peak flow meters. Also available for the Airzone.

Visit the ATM website or call today for complete information on our Asthma and COPD management products.

Distributed By:

Ordering Information

	Order #	Description/Quantity	NHRIC#
AllFLOW	3103105	Standard - each or 10/case	3 8462103105 0
Mini-Wright	3103104	Standard - each or 10/case	3 8462103104 0
	3104001	Low Range - each or 10/case	3 8462104001 1
Windmill Trainer	3104272	for use with AFS/Standard/ All Flow - each	
Airzone	3130045	each or 10/case	3 8462130045 0
Windmill Trainer	3104272	for use with Airzone - each	
Disposable Mouthpieces	3122200	Adult - 100/bag	3 8462122200 4
	3123200	Pediatric - 100/bag	3 8462123200 3
	3122064	One-way exp. - 200/bx	3 8462122064 2
Disposable Adapters	3122120	Pediatric/Small Mouth Adult - bag of 20	3 8462122120 5
	3122220	Pediatric/Small Mouth Adult - box of 100	3 8462122220 2



Alliance Tech Medical, Inc.
PO Box 6024 / Granbury, TX 76049
1.800.848.8923 / 817.326.6357
Fax: 817.326.2182 / www.alliancetekmedical.com

Eastern Office 302.838.2638 **Central Office** 800.848.8923
Southern Office 817.326.3183 **Western Office** 800.848.8923